



BUNNELL™ Splints INSTRUCTIONS #13-A & 13.1-A

COMBINATION OPPENHEIMER W/REVERSE KNUCKLE BENDER

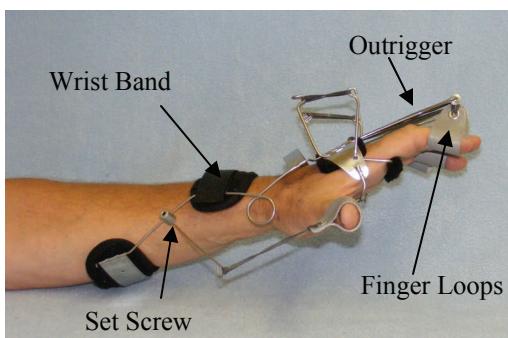

#13-A Coiled Spring Wire at Wrist

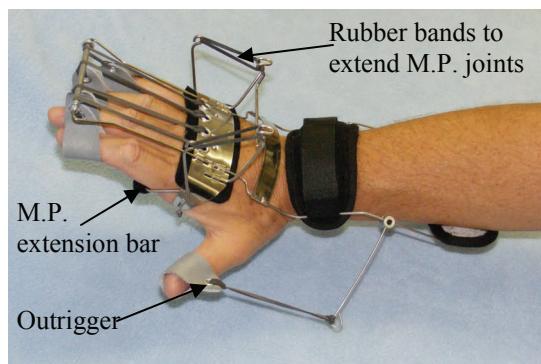

#13.1-A Straight Spring Wire at Wrist—Less extension force at wrist than #13-A

SET-UP AND ASSEMBLY

Position of the thumb outrigger determines if splint is for RIGHT or LEFT

To change position of the outrigger on the splint, loosen set screw that holds the outrigger to the lateral wire [hex key provided]. Reposition and tighten. The wire outrigger itself can be repositioned using the other set screw.


Be sure the wire end is flush with the aluminum housing when retightening.

Always remove splint from patient of all adjustments.


Attach finger loops to openings at top of hand band. Sizes medium, large, and extra large have four attachment openings (one for each finger). Extra small and small have two openings requiring two finger loops per opening.

ATTACH HAND AND FOREARM SECTIONS TOGETHER


Loosen the hand strap, reflect back, and expose openings for the forearm wire end on the forearm section. With the hand and forearm sections in proper position, insert the wire into the opening. Be sure the wire passes completely through the opening and rests in the channel.

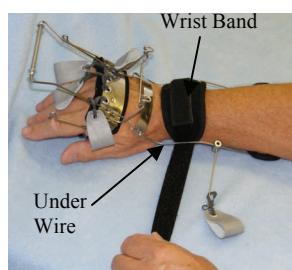

To close pass loop strap under lateral wire and back to secure with hook on top of wrist band. Wearing the wrist band/strap a little "looser" will apply less force to extend wrist.

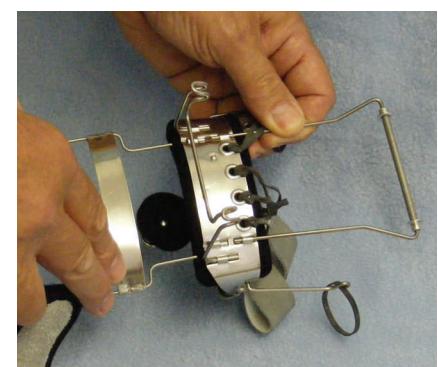

TO DON SPLINT Applying splint without rubber bands attached will facilitate proper positioning of splint. With splint on a flat surface open wrist band and slide hand into splint. Position M.P. extension bar under fingers between M.P. and P.I.P. joints. Position round pad in palm. Close wrist band. Attach rubber bands that extend M.P. joints to hooks. Place finger loops on fingers and thumb.

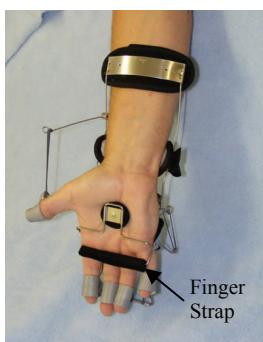

EVALUATE FIT, FUNCTION and INSTRUCT PATIENT ON USE
Splint is designed to extend wrist and finger joints while allowing active flexion.

TO ADJUST FORCE Remove splint and carefully bend lateral wires and/or change rubber bands.
Outrigger uses #31 rubber band

CLEAN with mild soap and damp towel. Remove liner, hand wash with cold water and mild detergent. Rinse thoroughly. Air dry completely before reapplication.

Avoid placing near open flame, heater, trunk of car, etc.

Immediately contact your health care professional if you are experiencing loss of sensation, swelling, skin irritations or have any questions or concerns.

CAUTION: Federal Law restricts this device to sale by or on the orders of a physician. To be dispensed by qualified healthcare professional.